

ZSI 375 PF ARTIFICIAL URINARY SPHINCTER

First Artificial Urinary Sphincter Pre-Filled



- Pre-connected / Pre-filled / 2 components device.
- One size fits all. Adjustable cuff from 3.75cm to 5cm (circumference).
- Ergonomic cuff moulded in Curve forme.
- Dissection-Insertion. No preparation.
- Reduced risk of infection due to short time in the OR.
- Possibility of increasing issued pressure.


Pre filled



Possibility of increasing
Issued pressure



1 3.75 to 5 cm adjustable
Cuff, moulded in a curve

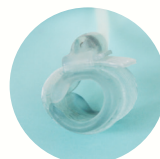

2 Deactivated when the button
(▼) portrudes at opposite
Side of the pump button



3 Pump button to empty
The cuff (press-release)




4 Activated when the button (▼)
portrudes at the side of the
Pump button



ZSI 475 HYDRAULIC PENILE IMPLANT




- All Standard sizes come in one box.
- Accessories for implantation are provided in the box / except dilators.
- Reduce time of kit preparation.
- PVP coating.
- All components are biocompatible.
- Larger diameter when inflated, resistive multiple folding, better rigidity.
- Extra large size 22cm available. *

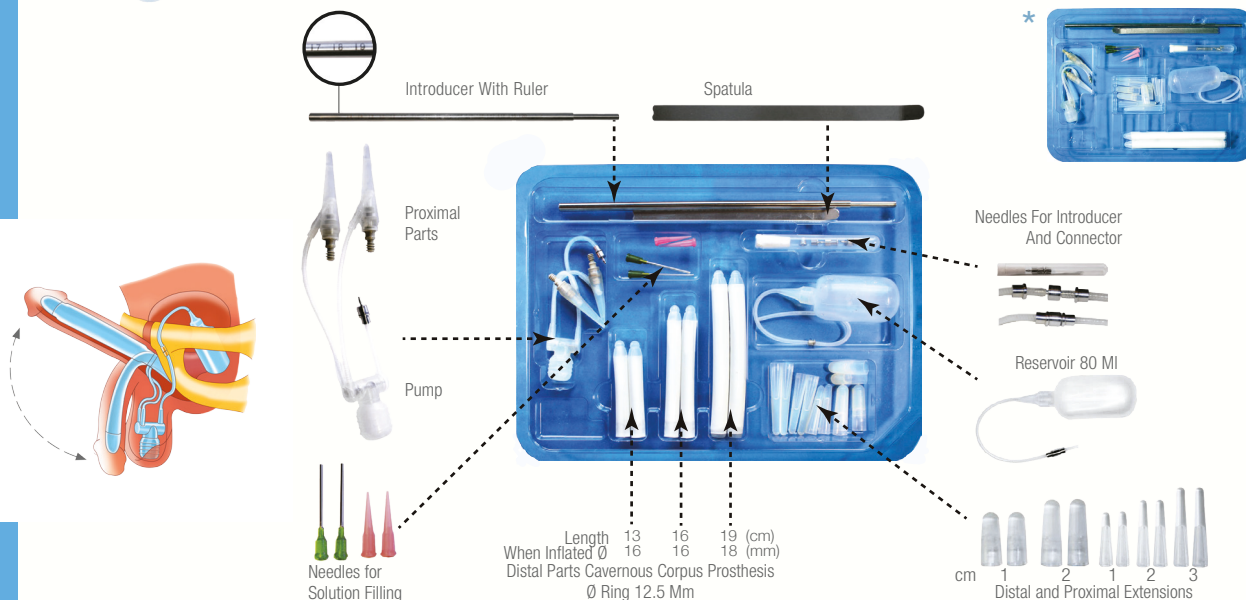


ZSI 100 MALLEABLE PENILE IMPLANT

- ZSI provide 2 bending rigidities.
Standard : designed for the best compromised.
Extra rigid : designed for strongest erection.
- Pair of implants fits all sizes. Adjustable from 12cm to 25cm.
- Adjustable distal and proximal for easy implantation.
- Biocompatible interior rod is made of Nitinol (Nickel-Titanium) for long term use.
- **GOOD ENDURANCE - GOOD RIGIDITY - GOOD PLASTICITY.**

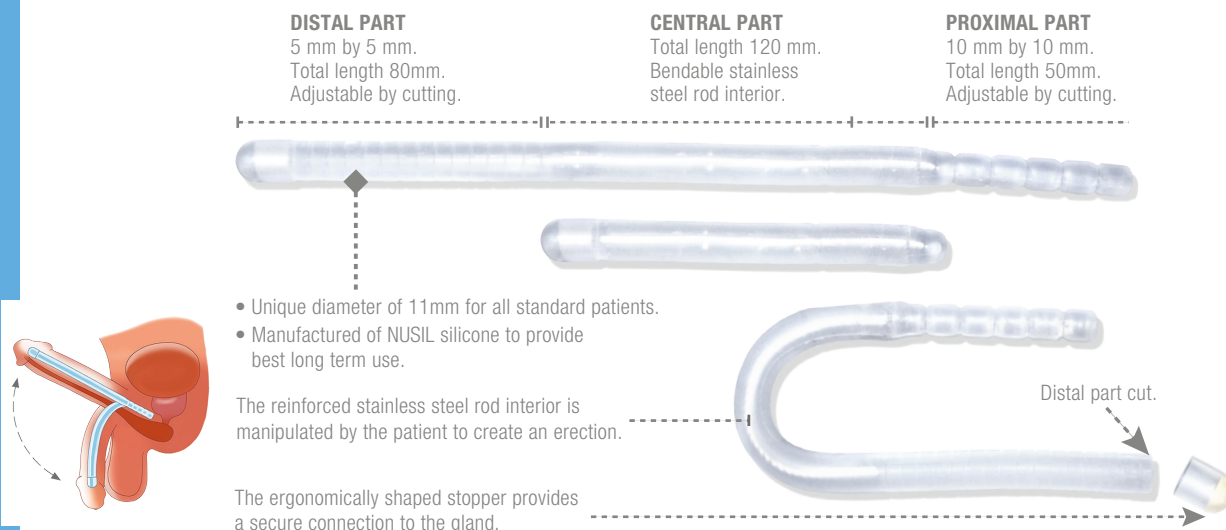

ZSI 100 FTM MALLEABLE PENILE IMPLANT FOR PHALLOPLASTY

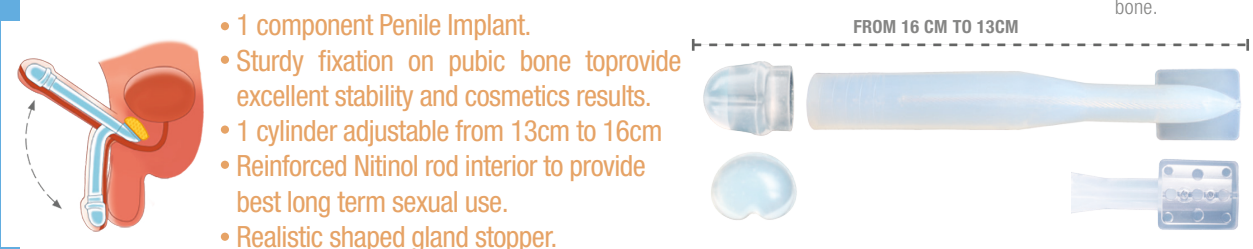

ZSI 475 FTM HYDRAULIC PENILE IMPLANT FOR PHALLOPLASTY

