


# ZSI 100 CF

## SOFT PENILE IMPLANT

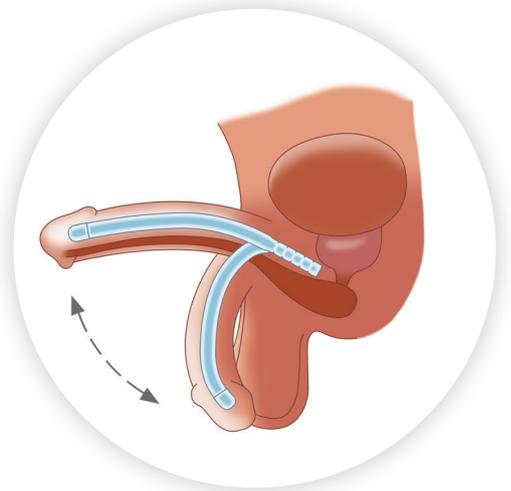

A soft solution to  
erectile dysfunction

- Pair of implants, one size fits all.
- For more information, visit [www.zsimplants.ch](http://www.zsimplants.ch)

# THE ZSI 100 CF SOFT PENILE IMPLANT DESIGNED TO PROVIDE EXCELLENT COSMETIC RESULTS

## DESCRIPTION:



The soft Penile Implant cable free ZSI 100 CF is designed to treat erectile dysfunction or to keep space free of fibrosis during healing after inflatable penile implant removal because of infection. It is composed of two flexible cylinders made of silicone with PVP coating to fight against infection. They are inserted in the corpus cavernosa.

## INSTRUCTION FOR USE

**1. ZSI 100 CF to treat erectile dysfunction:** the patient can get a flexible penis with erection stiff enough to penetrate the vagina.

**Procedure:** open the corpus cavernosa and make a space with dilator between inner corpus cavernosa wall and erectile tissue. Insert ZSI 100 CF into each corpus cavernous and close.

**2. ZSI 100 CF to keep a space into the corpus cavernous:** in case of inflatable penile implant removal for infection, the surgeon can keep the space into the corpus cavernous with the ZSI 100 CF. After 3 months or 6 months, the new inflatable penile implant insertion will be easier than to the ZSI 100 CF. Note that the ZSI 100 D9 CF/ ZSI 100 D11 CF/ ZSI 100 D13 CF get a PVP coating to fight against infection.

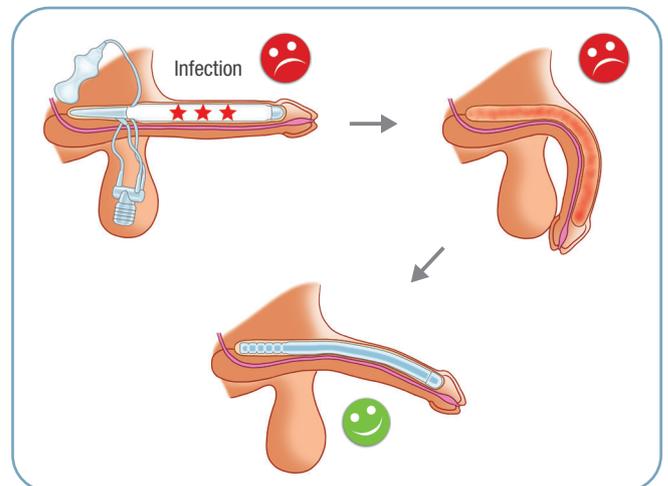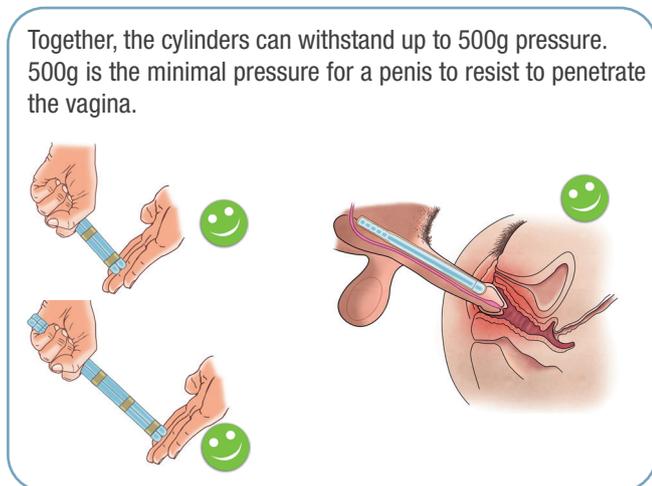

**PAIR OF IMPLANTS WITH PVP COATING, 3 DIAMETERS: 9 MM / 11MM / 13 MM**  
**ADJUSTABLE PROXIMALLY AND DISTALLY FROM 12 TO 25CM**  
**EASY IMPLANTATION AND OPTIMAL PATIENT COMFORT**

**PROXIMAL PART**  
 10 mm by 10 mm.  
 Total length 50mm.  
 Adjustable by cutting.

- 3 diameters: 9 mm / 11mm / 13 mm fit all standard patients.
- **NUSIL** silicone for optimal durability and biocompatibility.

